

***METHODS AND COMPOSITIONS FOR THE  
DETECTION OF MICROBIAL CONTAMINANTS***

Related Applications

This application claims priority to and the benefit of U.S.S.N. 60/455,632, filed March 17, 2003, the contents of which are incorporated by reference herein.

Field of the Invention

5           The present invention relates generally to methods and compositions for detecting and/or quantifying microbial contaminants in a sample. More particularly, the invention relates to methods and compositions using a hemocyte lysate for detecting and/or quantifying microbial contamination in a sample.

Background of the Invention

10           Microbial contamination by, for example, Gram positive bacteria, Gram negative bacteria, yeast, fungi, and molds may cause severe illness and, in some cases, even death in humans. Manufacturers in certain industries, for example, the pharmaceutical, medical device, and food industries, must meet exacting standards to verify that their products do not contain levels of microbial contaminants that would otherwise compromise the health  
15 of the recipient. These industries require frequent, accurate, and sensitive testing for the presence of such microbial contaminants to meet certain standards, for example, standards imposed by the United States Food and Drug Administration (USFDA) or Environmental Protection Agency. By way of example, the USFDA requires certain manufacturers of pharmaceuticals and invasive medical devices to establish that their  
20 products are free of detectable levels of Gram negative bacterial endotoxin.

Furthermore, when people become infected with Gram negative bacteria, the bacteria may produce and secrete fever-inducing bacterial endotoxins. Bacterial

endotoxins can be dangerous and even deadly to humans. Symptoms of infection may range from fever, in mild cases, to death. In order to promptly initiate proper medical treatment, it usually is important to identify, as early as possible, the presence of an endotoxin and, if possible, the concentration of the endotoxin in the patient.

5           To date, a variety of assays have been developed to detect the presence and/or amount of a microbial contaminant in a test sample. One family of assays use hemocyte lysates prepared from the hemolymph of crustaceans, for example, horseshoe crabs. These assays typically exploit, in one way or another, a clotting cascade that occurs when the hemocyte lysate is exposed to a microbial contaminant. A currently preferred  
10 hemocyte lysate is amebocyte lysate (AL) produced from the hemolymph of a horseshoe crab, for example, *Limulus polyphemus*, *Tachypleus gigas*, *Tachypleus tridentatus*, and *Carcinoscorpius rotundicauda*. Amebocyte lysates produced from the hemolymph of *Limulus*, *Tachypleus*, and *Carcinoscorpius* species are referred to as *Limulus amebocyte lysate* (LAL), *Tachypleus amebocyte lysate* (TAL), and *Carcinoscorpius amebocyte*  
15 *lysate* (CAL), respectively.

Routine assays that use LAL include, for example, gel clot assays, end point turbidometric assays, kinetic turbidometric assays, and endpoint chromogenic assays (Prior (1990) "Clinical Applications of the *Limulus* Amebocyte Lysate Test" CRC PRESS 28-34). These assays, however, suffer from one or more disadvantages including reagent  
20 expense, assay speed and limited sensitivity ranges. Also, these assays typically require that samples be sent to a testing facility removed from the origin of the sample being tested. As a result, it may take hours to weeks before a problem can be detected and remedied. Accordingly, there is an ongoing need for faster and more sensitive methods, and portable test systems employing such methods, that overcome the need to submit  
25 samples to an off-site testing facility.

### Summary of the Invention

The invention is based, in part, upon the discovery that it is possible to make and use optical cartridges containing an immobilized hemocyte lysate for use in hemocyte

lysate-based assays. These cartridges may be used alone or in combination with optical detectors, for example, hand held optical detectors, to permit the assay of samples in the field, thereby obviating the need to send samples to an off-site testing facility.

Accordingly, the cartridges can be used in a point-of-use test system. In addition, the  
5 invention is based, in part, upon the discovery of a rapid, sensitive, multi-step kinetic assay for determining the presence and/or amount of microbial contaminant in a sample of interest. This type of assay can be implemented in a cartridge, or in any other desirable assay format. In addition, the invention is based, in part, upon the discovery of a glucan-specific hemocyte lysate and a method of making such a lysate. The glucan-specific  
10 lysate may be incorporated into such a cartridge and/or used in a multi-step kinetic assay.

In one aspect, the invention provides a test cartridge for determining the presence and/or amount of a microbial contaminant in a sample. The cartridge comprises a housing defining at least one fluid inlet port, at least one optical cell, and at least one conduit having a fluid contacting surface that connects and thus provides fluid flow  
15 communication between the fluid inlet port and the optical cell. The cartridge further comprises a hemocyte lysate disposed on a first region of the fluid contacting surface of the conduit, so that when a sample is applied to the fluid inlet port, the sample traverses the region and solubilizes the hemocyte lysate during transport of the sample-lysate mixture to the optical cell. This type of cartridge can be used to perform, for example,  
20 endpoint turbidometric and kinetic turbidometric assays.

The cartridge optionally may further comprise a chromogenic substrate that acts as a substrate for one or more of the enzymes in the hemocyte lysate. The chromogenic substrate may be disposed in the first region, for example, combined with the hemocyte lysate. In this format, the sample resolubilizes or starts to resolubilize the hemocyte  
25 lysate and chromogenic substrate at substantially the same time. This type of cartridge can be used to perform, for example, a kinetic chromogenic assay. Alternatively, the chromogenic substrate may be disposed on a second region of the fluid contacting surface of the conduit, so that when the sample moves along the conduit toward the optical cell it contacts and reconstitutes the hemocyte lysate and chromogenic substrate at different

regions of the conduit. In one embodiment, the second region is located downstream of the first region (i.e., the second region is located closer to the optical cell than the first region). This type of cartridge can be used to perform, for example, endpoint chromogenic and multi-step kinetic chromogenic assays, as discussed in more detail  
5 herein below.

In addition, a pre-selected amount of an agent representative of a microbial contaminant, or “spike”, such as a bacterial endotoxin, a (1→3)-*B*-D glucan, or other microbial cell wall constituent, is disposed on a region of the fluid contacting surface of the conduit. The inclusion of the agent or spike is particularly useful as it provides a  
10 positive control to demonstrate that an assay is working, and can also demonstrate whether an inhibitor or enhancer is present in the sample. The agent or spike may be disposed on the first region, the second region, or on another region of the conduit.

It is contemplated that the number of fluid inlet ports, optical cells, and conduits in a particular cartridge may vary depending on the number of samples or microbial  
15 contaminants being tested at a particular time. A cartridge may be used to simultaneously assay duplicates of the same sample or simultaneously assay two or more different samples of interest. Alternatively, two or more different hemocyte lysates may be disposed on the cartridge so that it is possible to determine the presence and/or amount of two or more different microbial contaminants at substantially the same time. In addition,  
20 several chromogenic substrates with different enzyme specificities and optical properties, for example, light absorption transmission, and/or fluorescent properties, may be applied to the cartridge. This may permit the detection of two or more different microbial contaminants at substantially the same time.

In one embodiment, the cartridge comprises a housing that defines (i) a first fluid  
25 inlet port, a first optical cell, and a first conduit having a fluid contacting surface that connects and thus provides fluid flow communication between the first fluid inlet port and the first optical cell, and (ii) a second fluid inlet port, a second optical cell, and a second conduit having a fluid contacting surface that connects and thus provides fluid flow communication between the second fluid inlet port and the second optical cell.

Within the housing, a first hemocyte lysate is disposed on a first region of the fluid contacting surface of the first conduit, so that when a sample is applied to the first fluid inlet port, the sample traverses the region and reconstitutes and/or solubilizes the first hemocyte lysate during transport to the first optical cell. Also within the housing, a  
5 second hemocyte lysate is disposed on a first region of the fluid contacting surface of the second conduit, so that when a sample is applied to the second fluid inlet port, the sample traverses the region and reconstitutes and/or solubilizes the second hemocyte lysate during transport to the second optical cell.

In one embodiment, a chromogenic substrate is disposed on a second region of the  
10 fluid contacting surface of the first conduit and/or a chromogenic substrate is disposed on a second region of the fluid contacting surface of the second conduit. In each embodiment, the second region preferably is located downstream of the first region in each conduit. In another embodiment, different chromogenic substrates may be disposed on fluid contacting surfaces of the first and second conduits so that two different reactions  
15 may be monitored in the same cartridge.

In another embodiment, a pre-selected amount of an agent representative of a microbial contaminant or spike is disposed on the fluid contacting surface of the first or second conduit. The spike may be disposed on the first region or on another region of the conduit. The spike may be useful as a positive control for the assay (i.e., indicates  
20 whether a valid test was run), and may also provide information on whether an enhancer or inhibitor is present in the sample.

As will be discussed in more detail below, the cartridges of the invention may be adapted for use in a variety of different assays. The presence of a microbial contaminant may be indicative of, for example, past or present bacterial, yeast, fungal, or mold  
25 infection. It is contemplated that, by using the appropriate hemocyte lysate, the cartridge may be used to detect the presence of a bacterial, yeast, fungal, or mold contaminant in a sample. The cartridges of the invention are particularly useful at detecting the presence, and/or determining the amount, of a Gram negative bacterial endotoxin or glucan in a sample.

During use, a sample to be tested is introduced into the sample inlet port of the cartridge and is allowed to move to the optical cell. Movement of the sample along the conduit can be passive or can be induced by an external force, for example, via positive or negative pressure in the conduit. For example, the sample can be pulled along the  
5 conduit to the optical cell via suction induced by a pump connected to a pump port located downstream of the optical cell. A change in an optical property of the sample is detected in the optical cell, the change being indicative of the presence of a microbial contaminant in the sample. In addition, the time in which a pre-selected change occurs in an optical property of the sample can be determined and compared against a  
10 predetermined standard curve to determine the concentration of the microbial contaminant in the sample. The optical property may include a color change, a change in absorbance or transmittance, a change in turbidity, a change in fluorescence or other change that can be detected in a detector, spectrophotometer or the like. The change in the optical property may be, for example, an increase in absorbance of light of a pre-  
15 selected wavelength, or may be a decrease in transmission of light of a pre-selected wavelength. As discussed below, the cartridge may be adapted to perform any one of a number of endpoint or kinetic chromogenic, or turbidometric assays.

In another aspect, the invention provides methods of preparing the cartridge by disposing, for example, by drying, a hemocyte lysate onto a solid surface of at least one  
20 conduit of the cartridge. The hemocyte lysate may then be reconstituted into an active form upon resolubilization of the hemocyte lysate. A volume of a mixture comprising a hemocyte lysate of interest and a resolubilizing agent and/or an anti-flaking agent is applied to the surface of at least one conduit and dried. The hemocyte lysate used will depend upon the assay for which the cartridge will be used. The resolubilizing agent is an  
25 agent that, either alone or in combination with another resolubilizing agent, facilitates the resolubilization of one or more components of the hemocyte lysate once the lysate is exposed to a fluid sample. The resolubilizing agent preferably also stabilizes the lysate in its dried form. The resolubilizing agent provided in the mixture facilitates the stability of the reagents and their dissolution during the assay. Resolubilizing agents include, for  
30 example, one or more sugars, salts, or combinations thereof. Preferred sugar

resolubilizing agents include, for example, mannitol, mannose, sorbitol, trehalose, maltose, dextrose, sucrose, and other monosaccharides or disaccharides. The anti-flaking agent included in the mixture further stabilizes the reagents and reduces flaking of the dried lysate. The anti-flaking agent preferably also stabilizes the lysate in its dried form.

5 Preferred anti-flaking agents include, for example, one or more polymers, for example, polyethylene glycol, polyvinyl pyrrolidone, polyvinyl alcohol, mannitol, dextran, and proteins, for example, serum albumin.

In one embodiment, the lysate/resolubilizing agent/anti-flaking agent mixture is disposed in a first region of at least one conduit of the cartridge. The mixture then is  
10 dried onto a surface of the conduit in an environment having a temperature of about 4°C to about 40°C, more preferably, from about 10°C to about 35°C, more preferably, from about 15°C to about 30°C and a relative humidity of about 0% to about 30%, more preferably, from about 2% to about 20%, more preferably, from about 4% to about 10%. Preferably, the temperature is about 25°C and the relative humidity is about 5%. Drying  
15 preferably occurs for about 10 minutes to about 8 hours, more preferably for about 1 hour in a temperature regulated drying chamber. The lysate/resolubilizing agent/anti-flaking agent mixture optionally may also include an agent representative of a microbial contamination (i.e., a spike), for example, a bacterial endotoxin, a (1→3)-*B*-D glucan or other microbial cell wall constituents.

20 In another embodiment, the mixture is dried onto the surface of the conduit by lyophilization or freeze drying, for example, at temperatures below 0°C, for example, from about -75°C to about -10°C, more preferably from about -30°C to about -20°C.

In addition, a chromogenic substrate, comprising a resolubilizing agent and an anti-flaking agent is applied to a second region of the cartridge and dried onto the  
25 cartridge as described above. The chromogenic substrate may comprise an anti-frothing agent, for example, polyvinyl alcohol and polypropylene glycol.

Although the drying procedure is discussed in relation to a test cartridge, the drying procedure can be used to dry the lysate onto a variety of different solid supports.

For example, the method can also be used to dry the lysate on the surface of a well disposed within or defined by a solid support, for example, a 12-well or a 96-well plate.

In another aspect, the invention provides an improved method referred to as a two-step or multi-step kinetic assay for detecting the presence and/or quantifying the amount of a particular microbial contaminant in a sample. The sample to be tested is contacted with a hemocyte lysate comprising an activatable enzyme, for example, a pro-clotting enzyme or a clotting enzyme, that is activated if the microbial contaminant is present in the sample. After contacting the sample with the lysate, the sample-lysate mixture is incubated for a preselected period of time. Then, the sample-lysate mixture is contacted with a substrate, for example, a chromogenic substrate, for the activated enzyme. If the sample-lysate substrate mixture contains an activated enzyme, the activated enzyme produces a change in the substrate, which in turn produces a predetermined change in an optical property of the sample-lysate-substrate mixture. The time in which the predetermined change occurs then is determined. The resulting time then is compared to a predetermined standard curve to determine whether the microbial contaminant is present in the sample and/or to determine the amount of microbial contaminant present in the sample. For example, the concentration of microbial contaminant in a sample can be measured by comparing the time required to produce the predetermined change in the optical property against a predetermined standard curve of the microbial contaminant. Using this type of assay, it is possible to adjust the timing of the steps to produce an assay of predetermined sensitivity and duration.

The optical property measured can be a change (e.g., an increase or decrease) in an optical property, for example, absorbance at a particular wavelength, transmittance at a particular wavelength, fluorescence at a particular wavelength, or optical density. For example, the optical property may be a change in absorbance or transmittance at a wavelength in the range from about 200 nm to about 700 nm, and more preferably in the range from about 350 nm to about 450 nm.

The chromogenic substrate may be any substrate for a lysate enzyme that is activated (e.g., hydrolyzed) to cause a detectable chromogenic or fluorogenic change, for



example, by release of a chromophore or a fluorophore, that is detectable by an optical detector. In one embodiment, the chromogenic substrate for LAL contains a para-nitroaniline chromophore, such as that, for example, in the chromogenic substrate acetate-Ile-Glu-Ala-Arg-pNA. The proteases in the LAL cleave colorless tetrapeptide to release  
5 the pNA group, which causes a color change. Cleavage of the tetrapeptide simulates the cleavage reaction of the proteases in the LAL with coagulogen, a clotting component that contains the tetrapeptide. As a result, by using this chromophore, it is possible to measure the progress of the reaction by measuring the change in optical density at about 395 nm. Other chromophores may include dinitrophenyl alanine, cyclohexyl alanine and  
10 the like. Alternatively, the substrate may contain a fluorophore, for example, 7-amino-4-methyl coumarin, 7-amino-4-trifluoromethyl coumarin, and 4-methoxy-2-naphthylamine. Fluorogenic substrates for LAL that contain N-methylcoumarin as a leaving group are available from Enzyme Systems Products, Livermore, CA.

The multi-step kinetic assay may be performed in a variety of formats, for  
15 example, in a tube, cuvette, cartridge, well on a solid support (such as a 96-well multi-well plate), or other vessel suitable for use in combination with an optical detector, for example, a spectrophotometer, fluorimeter, luminometer, or the like.

In another aspect, the invention provides a method for producing an amebocyte lysate depleted of Factor C activity. The method comprises the steps of: (a) providing a  
20 preparation of amebocytes; and (b) lysing the amebocytes in the presence of at least 0.05M salt to provide an amebocyte lysate preparation depleted of Factor C activity. The method optionally includes the step of, after step (b), removing cellular debris, for example, cell membranes, and then harvesting the remaining lysate. The cellular debris may be sedimented by centrifugation and the remaining supernatant harvested.

25 In one embodiment, the salt may comprise a monovalent cation, for example, a sodium or potassium salt. Salts useful in the practice of the invention include, for example, sodium chloride, potassium chloride, sodium acetate, and potassium acetate. The salt concentration can be in the range from 0.15 M to about 6 M, more preferably from about 0.25 M to about 4 M, and more preferably from about 1 M to 2 M. However,

the precise concentration of salt necessary to remove or reduce Factor C activity may be determined by routine experimentation. For example, amebocytes can be lysed in the presence of difference concentrations of salt, and the residual lysates can then be checked to see whether a coagulin clot forms in the presence of a bacterial endotoxin. The  
5 foregoing method can produce a glucan-specific amebocyte lysate that is substantially free of Factor C activity. The lysate, therefore, retains Factor G activity but is depleted of Factor C activity.

In another aspect, the invention provides an amebocyte lysate substantially free of Factor C activity, wherein the lysate comprises at least about 0.25M salt and wherein the  
10 lysate is capable of producing a coagulin gel in the presence of glucan. The lysate may comprise, from about 0.25 M salt to about 6 M salt, from about 0.5 M salt to 4 M salt, and from about 1 M to about 2 M salt. In one embodiment, the salt contains a monovalent cation, for example, a sodium ion or a potassium ion. For example, the salt may include sodium chloride, potassium chloride, sodium acetate, potassium acetate, or a  
15 combination thereof.

The foregoing and other objects, features and advantages of the present invention will be made more apparent from the following drawings and detailed description of preferred embodiments of the invention.

#### Brief Description of the Drawings

20 The objects and features of the invention may be better understood by reference to the drawings described below in which,

Figure 1 is a schematic representation of the coagulation system present in amebocytes;

25 Figures 2A-2D are schematic illustrations of an exemplary cartridge of the invention in perspective view (Figure 2A), top view (Figure 2B), side view (Figure 2C), and end view (Figure 2D);

Figure **3A-3B** are schematic illustrations of an exemplary cartridge of the invention wherein each conduit has a separate fluid inlet port (Figure **3A**), and wherein two conduits share a single common fluid inlet port (Figure **3B**);

5 Figures **4A-4B** are schematic illustrations of an exemplary cartridge of the invention wherein each conduit has its own fluid inlet port (Figure **4A**), and wherein a pair of conduits share a single common fluid inlet port (Figure **4B**);

10 Figures **5A-5D** are schematic illustrations of an exemplary cartridge in which Figure **5A** is a view of a bottom half of an exemplary cartridge of the invention showing the locations of immobilized hemocyte lysate and chromogenic substrate, Figure **5B** is a view of a top half of an exemplary cartridge of the invention showing the location of an immobilized agent representative of a microbial contaminant (i.e., spike), Figure **5C** is a cross-sectional view of the fabricated cartridge through section A-A', and Figure **5D** is a cross-sectional view of the fabricated cartridge through section B-B';

15 Figure **6A-6D** are schematic illustrations of two exemplary cartridges, in which Figure **6A** is a top view of a first embodiment of a cartridge, Figure **6B** is a side view of the first embodiment of the cartridge pictured in Figure **6A**, Figure **6C** is a top view of a second, different embodiment of a cartridge, and Figure **6D** is a side view of the second embodiment of the cartridge of Figure **6C**;

20 Figure **7** is a flow chart for an exemplary multi-step kinetic chromogenic assay of the invention;

Figure **8A-8B** are schematic illustrations showing a cartridge in combination with an exemplary optical detector in which Figure **8A** shows the cartridge being inserted into the detector, and Figure **8B** shows the cartridge actually inserted into the detector;

25 Figure **9** is a graphical representation showing changes in light absorbance or optical density (full line) or light transmittance (dashed line) through an optical cell of an exemplary cartridge of the invention during a multi-step kinetic chromogenic assay;

Figures **10A-10B** are graphical representations of absorbance values in an endpoint chromogenic assay where Figure **10A** shows the absorbance values of endotoxin standards (1.0 Endotoxin Units (EU)/mL (A1), 0.5 EU/mL (A2), 0.25 EU/mL (A3), and 0.125 EU/mL (A4)) over time in an endpoint chromogenic assay performed in a cartridge of the invention, and Figure **10B** shows a standard curve generated by plotting the absorbance values of each concentration of endotoxin at T= 780 seconds;

Figures **11A-11D** are graphical representations of absorbance values for a kinetic chromogenic assay performed in a cartridge of the invention, where Figure **11A** shows the absorbance values of a 5.0 EU/mL endotoxin standard, Figure **11B** shows the absorbance values of a 0.5 EU/mL endotoxin standard, Figure **11C** shows the absorbance values of a 0.05 EU/mL endotoxin standard, and Figure **11D** is a standard curve generated by plotting the log of endotoxin concentration (X-axis) versus the log of absorbance value at onset time (Y axis);

Figure **12** is a graphical representation of a standard curve for a multi-step kinetic chromogenic assay performed in a cartridge of the invention, generated by plotting the log of endotoxin concentration (X-axis) versus the log of absorbance value at onset time (Y axis);

Figure **13** is a graphical representation of a standard curve for a single-step kinetic chromogenic assay performed in a microtiter plate, generated by plotting the log of endotoxin concentration (X-axis) versus the log of absorbance value at onset time (Y axis);

Figures **14** is a graphical representation of a standard curve for an endpoint chromogenic assay performed in a microtiter plate, generated by plotting the absorbance values (Y axis) of each concentration of endotoxin (X axis); and

Figure **15** is a graphical representation of a standard curve for a multi-step kinetic chromogenic assay performed in a microtiter plate, generated by plotting the log of

endotoxin concentration (X-axis) versus the log of absorbance value at onset time (Y axis).

Figures 16A-B are graphical representations of amebocyte lysate kinetic reactions for standard LAL and glucan-specific LAL. Figure 16A represents assays performed using lipopolysaccharide, and Figure 16B represents assays performed using glucan. Row A in each Figure represents standard LAL containing both the Factor C and Facion G cascades, Row B in each Figure represents glucan-specific lysate (i.e., Factor G-specific lysate) prepared by lysing amebocytes in 1M sodium chloride and Row C represents glucan-specific lysate prepared by lysing amebocytes in 2M NaCl. Columns 1-5 represent several dilutions of the lipopolysaccharide (10 ng/ml, 1ng/ml, 100 pg/ml, 10 pg/ml, 0) or glucan (100 µg/ml, 10 µg/ml, 1 µg/ml, 100 ng/ml, 0) added to each sample, wherein the concentration decreases from column 1 to column 5;

Figure 17 is a graphical representation of a logarithmic plot of a glucan standard curve obtained using a multi-step kinetic assay;

Figure 18 is a graphical representation of an assay measuring the fluorescence emitted from different concentrations of a fluorogenic substrate Glu-Gly-Arg-AMC in a multi-step kinetic assay in the cartridge format;

Figure 19 is a graphical representation showing the proportional displacement of various concentrations of lipopolysaccharide-fluorescein isothiocyanate (1 µg/ml, 100 ng/ml, 10 ng/ml, 1 ng/ml, 100 pg/ml, 10 pg/ml, and control) by different dilutions of a fluorescein-labeled ligand (1/50, 1/150, and 1/450).

In the drawings, which are not necessarily drawn to scale, like characters refer to the same or similar parts throughout the Figures.

#### Detailed Description of the Invention

The invention provides an optical cartridge containing an immobilized hemocyte lysate for use in hemocyte-lysate based assays. These cartridges may be used alone or

together with an optical detector, for example, a hand held optical detector. In addition, the invention provides a rapid, sensitive, broad range, multi-step assay that is useful in determining the presence and/or amount of a microbial contaminant in a sample.

Although the cartridge and method may be used separately, they are particularly effective  
5 when combined together to provide a system that can be used in the field to provide rapid test results. This facilitates quicker elimination and/or treatment of microbial contamination. In addition, the invention provides a Factor C specific lysate for detecting the presence and/or amount of glucan in a sample. The lysate, therefore, can be used to determine the presence and/or amount of a yeast or mold contaminant in a sample.

#### 10 The Cartridge

It is contemplated that the cartridges of the invention may be formulated with one or more hemocyte lysates and used in a variety of assays to detect the presence and/or amount of a microbial contaminant in a sample. A number of hemocyte lysate-based assays for the detection and/or quantification of a microbial contaminant can be  
15 performed in the cartridge of the invention, for example, as illustrated in Figure 2. The cartridge may be used on its own and the test result detected by eye or may be used in combination with an optical detector, for example, a hand-held optical detector as shown and described in U.S. Patent No. Des. 390,661.

By way of example and as illustrated in Figures 2A-2D, cartridge 1 has a  
20 substantially planar housing fabricated, for example, from a moldable biocompatible material. The housing may be fabricated from any material, however, transparent and/or translucent glass or polymers are preferred. Preferred polymers include, for example, polystyrene, polycarbonate, acrylic, polyester, optical grade polymers, or any plastic such that the optical cell is substantially transparent. The housing contains at least one fluid  
25 inlet port 4, at least one optical cell 6, and at least one conduit 8 having a fluid contacting surface for providing fluid flow communication between the fluid inlet port 4 and optical cell 6. The only requirements for the optical cell 6 are that it defines a void capable of containing a sample to be tested and that a portion of the optical cell 6 is transparent to light. Cartridge 1 may also have at least one pump port 12 in fluid flow communication

with fluid inlet port 4 and optical cell 6 for attaching the cartridge 1 to a pump. The pump may then impart a negative pressure via pump port 12 to pull the sample from fluid inlet port 4 to optical cell 6. A hemocyte lysate is disposed on a first region 14 of the fluid contacting surface of conduit 8, so that when a sample is applied to fluid inlet port 4, the sample traverses region 14 and solubilizes or reconstitutes the hemocyte lysate into the sample as it moves toward optical cell 6. This type of cartridge 1 may be used for performing, for example, an endpoint turbidometric or a kinetic turbidometric assay. In an embodiment, a chromogenic substrate may also optionally be applied to the surface of the conduit 8 at first region 14 together with the hemocyte lysate. This type of cartridge 1 may be used for performing, for example, a kinetic chromogenic assay.

In a preferred embodiment, as illustrated in Figures 2A-2D, a second region 16 of the fluid contacting surface of conduit 8 is spaced apart from and downstream of first region 14. In this configuration, hemocyte lysate is disposed at first region 14 and a chromogenic substrate is disposed at second region 16, so that after the sample is contacted with the hemocyte lysate in region 14, the sample-lysate mixture traverses conduit 8 and contacts the chromogenic substrate in region 16. The sample-lysate-substrate mixture then traverses conduit 8 to optical cell 6. This type of cartridge may be used for performing, for example, an endpoint chromogenic assay or a multi-step kinetic chromogenic assay, as discussed in more detail below.

Depending upon the type of assay to be performed, a pre-selected amount of an agent representative of a microbial contaminant, or "spike," such as a bacterial endotoxin, is disposed on first region 14 of the fluid contacting surface of one or more conduits 8. Alternatively, the spike may be disposed on a different region of the conduit 8.

It is contemplated that the cartridge 1 may have a variety of different configurations, demonstrated, for example, in Figures 3 and 4. Figure 3A shows a cartridge 1 comprising two conduits 8 with each conduit 8 having its own sample inlet port 4. Figure 3B shows a cartridge 1 comprising two conduits 8 with each conduit 8 sharing a common sample inlet port 4. Figure 4A shows a cartridge 1 comprising four

separate conduits 8 with each conduit 8 having its own sample inlet port 4. Figure 4B shows a cartridge 1 comprising two pairs of conduits 8, with each pair having its own sample inlet port 4.

The cartridges can be designed and used according to the type and/or number of tests required. For example, a single sample may be tested, for example, in duplicate or triplicate, for example, for research laboratory uses or for medical device and biopharmaceutical testing. Alternatively, two or more different samples may be tested individually, for example, for dialysis facility testing of water and dialysate. The cartridge preferably is a single-use, disposable cartridge that is discarded after one use.

10 The cartridge of the invention can use approximately 20-100 fold less hemocyte lysate per sample than is used in the conventional endpoint chromogenic or kinetic chromogenic assays performed in multi-well plates, and thus provides a less costly and environmentally-friendlier test.

Once a particular assay format has been chosen, the cartridge may be fabricated as discussed below.

#### Cartridge Fabrication

All the reagents and materials used to prepare the cartridge preferably are free of the microbial contaminant for which the cartridge ultimately will be used to test.

It is contemplated that the cartridge may be fabricated with any hemocyte lysate of interest. As used herein, the term, "hemocyte lysate" is understood to mean any lysate or a fraction or component thereof, produced by the lysis and/or membrane permeabilization of hemocytes, for example, amebocytes and hemolymph cells, (i) extracted from a crustacean or insect and/or (ii) cultured *in vitro* after extraction from the host. Hemocyte cellular material that has been extruded from hemolymph cells by contact with a

20 membrane permeabilization agent such as a Ca<sup>2+</sup> ionophore or the like (i.e., extruded other than by lysis) or otherwise extracted without cellular lysis is also considered to be a

25 hemocyte lysate. A preferred hemocyte lysate is an amebocyte lysate prepared from the



blood of a crustacean, for example, a horseshoe crab or Jonah crab. It is also contemplated that hemocyte lysate may include a cocktail of one or more natural (e.g., purified) or synthetic components of the enzyme cascades shown in Figure 1.

As used herein, the term “amebocyte lysate” is understood to mean any lysate or  
5 fraction or component thereof produced by the lysis, extrusion, or extraction of the cellular contents from amebocytes extracted from a crustacean, for example, a horseshoe crab. The amebocyte lysate comprises at least one component of an enzymatic cascade (for example, as shown in Figure 1) and/or produces a clot in the presence of an endotoxin, for example, a Gram negative bacterial endotoxin and/or a glucan, for  
10 example, a (1→3)-β-D glucan, produced by a yeast or a mold. Preferred amebocyte lysates can be derived from horseshoe crabs, which include crabs belonging to the *Limulus* genus, for example, *Limulus polyphemus*, the *Tachypleus* genus, for example, *Tachypleus gigas*, and *Tachypleus tridentatus*, and the *Carcinoscorpius* genus, for example, *Carcinoscorpius rotundicauda*.

15 Crude lysates may be produced using the procedure as originally described in Levin *et al.* (1968) THROMB. DIATH. HAEMORRH. 19: 186, with modification, or in Prior (1990) “Clinical Applications of the *Limulus* Amebocyte Lysate Test” CRC PRESS 28-36 and 159-166, and in U.S. Patent No. 4,322,217. Other lysates may include those, for example, described in U.S. Patent Nos. 6,270,982 and 6,391,570.

20 Presently, LAL is employed as the amebocyte lysate of choice in many bacterial endotoxin assays because of its sensitivity, specificity, and relative ease for avoiding interference by other components that may be present in a sample. LAL, when combined with a sample containing bacterial endotoxin and optionally with certain LAL substrates, reacts with the endotoxin in the sample to produce a detectable product, such as a gel,  
25 increase in turbidity, or a colored or light-emitting product, in the case of a synthetic chromogenic substrate. The product may be detected, for example, either visually or by the use of an optical detector.

As shown in Figure 1, the coagulation system of hemolymph, like the mammalian blood coagulation system, comprises at least two coagulation cascades that include an

endotoxin-mediated pathway (the Factor C pathway) and a (1→3)-*B*-D glucan-mediated pathway (the Factor G pathway). See, for example, Morita *et al.* (1981) FEBS LETT. 129: 318-321 and Iwanaga *et al.* (1986) J. PROTEIN CHEM. 5: 255-268.

The endotoxin-mediated activation of LAL is well understood and has been  
5 thoroughly documented in the art. See, for example, Levin *et al.* (1968) *supra*; Nakamura  
*et al.* (1986) EUR. J. BIOCHEM. 154: 511; Muta *et al.* (1987) J. BIOCHEM. 101: 1321; and  
Ho *et al.* (1993) BIOCHEM. & MOL. BIOL. INT. 29: 687. When bacterial endotoxin is  
contacted with LAL, the endotoxin initiates a series of enzymatic reactions, referred to in  
the art as the Factor C pathway, that can involve three serine protease zymogens called  
10 Factor C, Factor B, and pro-clotting enzyme (see, Figure 1). Briefly, upon exposure to  
endotoxin, the endotoxin-sensitive factor, Factor C, is activated. Activated Factor C  
thereafter hydrolyses and activates Factor B, whereupon activated Factor B activates  
proclotting enzyme to produce clotting enzyme. The clotting enzyme thereafter  
hydrolyzes specific sites, for example, Arg<sup>18</sup>-Thr<sup>19</sup> and Arg<sup>46</sup>-Gly<sup>47</sup> of coagulogen, an  
15 invertebrate, fibrinogen-like clottable protein, to produce a coagulin gel. See, for  
example, U.S. Patent No. 5,605,806.


It is possible to produce an endotoxin-specific lysate by removing Factor G  
activity from the lysate, such as the Factor G depleted lysates produced by the methods  
described in U.S. Patent Nos. 6,391,570 and 6,270,982. Also, it is contemplated that  
20 lysates may be depleted of Factor G activity by the addition of certain inhibitors or  
modulators of Factor G activity, for example, certain detergents, saccharides,  
polysaccharides, and other reagents described in U.S. Patent Nos. 5,155,032; 5,179,006;  
5,318,893; 5,474,984; 5,641,643; 6,270,982; and 6,341,570. An endotoxin-specific  
lysate is a lysate capable of reacting with a bacterial endotoxin but in which the reactivity  
25 to (1→3)-*B*-D glucan has been depleted by at least 80%, more preferably by at least 90%,  
and more preferably by at least 95% relative to the crude lysate from which the  
endotoxin-specific lysate was prepared.

(1→3)-*B*-D glucans and other LAL reactive glucans, produced by microorganisms  
such as yeasts and molds, can also activate the clotting cascade of LAL, through a

different enzymatic pathway, referred to in the art as the Factor G pathway (see, Figure 1). Factor G is a serine protease zymogen that becomes activated by (1→3)- β-D glucan or other LAL reactive glucans. Upon exposure to (1→3)- β-D glucan, for example, Factor G is activated to produce activated Factor G. Activated Factor G thereafter  
5 converts the proclotting enzyme into clotting enzyme, whereupon the clotting enzyme converts coagulogen into coagulin.

As used herein, the term, “(1→3)-β-D glucan” is understood to mean any water soluble polysaccharide, disaccharide or derivative thereof that is (i) capable of inducing formation of a coagulin clot in crude *Limulus* amebocyte lysate, and (ii) contains at least  
10 two β-D glucosides, connected by a (1→3)-β-D glycosidic linkage (see Formula I). It is contemplated that such a polysaccharide or derivative thereof, in addition to containing a (1→3)-β-D glycosidic linkage, may also contain glucoside moieties connected by a variety of other glycosidic linkages, for example, via a (1→4)-β-D glycosidic linkage and/or by a (1→6)-β-D glycosidic linkage. It is contemplated that such (1→3)-β-D  
15 glucans may be isolated from a variety of sources including, without limitation, plants, bacteria, yeast, algae, and fungi, or alternatively may be synthesized using conventional sugar chemistries.

**Formula I**

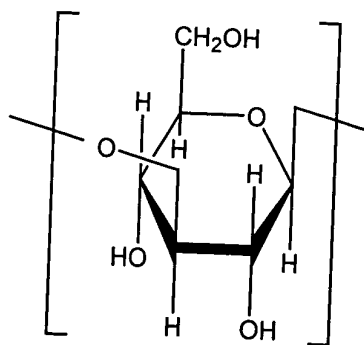

20

It is possible to produce a (1→3)-β-D glucan specific lysate by producing a lysate depleted of Factor C activity. As shown herein, it is possible to produce a glucan-specific

lysate by lysing amebocytes in the presence of at least 0.15 M salt, more preferably 0.25 M salt, for example, a salt containing a monovalent cation, such as sodium or potassium ions. For example, the amebocytes are lysed in about 0.15 M to about 6 M salt, for example, sodium chloride. Alternatively, the amebocytes are lysed in a solution  
5 containing from about 0.25 M to about 4 M salt, or from about 1 M to about 2M salt. When using sodium chloride, it appears that the amoebocyte preparation loses substantial Factor C activity when the amebocytes are lysed in a solution containing 0.25 M sodium chloride. However, the concentration of other salts necessary to produce a comparable results may be determined by routine titration experiments. For example, amebocytes  
10 may be lysed in different concentrations of salt and the resulting lysates examined for their ability to produce a coagulin gel in the presence of a Gram negative bacterial endotoxin. The concentration of salt may be chosen where the resulting lysate has lost a substantial amount of reactivity to the bacterial endotoxin. Other salts that may be used include, but are not limited to, monovalent ionic salts, such as, potassium chloride,  
15 potassium acetate and sodium acetate. An exemplary method for producing a glucan specific lysate is described in Example 4. A glucan-specific lysate is a lysate capable of reacting with glycan, for example, (1→3)-β-D glucan, but in which reactivity to a bacterial endotoxin or lipopolysaccharide has been depleted by at least 80%, more preferably at least 90%, and more preferably at least 95% relative to the crude lysate from  
20 which the glucan-specific lysate was prepared.

Methods for enhancing the sensitivity of hemocyte lysate for endotoxin, for example, may include, without limitation, aging the crude hemocyte lysate, adjusting pH, adjusting the concentration of divalent cations, adjusting the concentration of coagulogen, chloroform extraction, and the addition of serum albumin, biocompatible buffers and/or  
25 biological detergents.

As will be apparent to one of ordinary skill, divalent metal salts, which are known to promote activation of the pro-clotting enzyme of hemocyte lysate, as well as buffers to avoid extremes of pH that could inactivate the clotting enzyme preferably are included in the lysate. Any of the buffers and salts that are understood in the art to be compatible

with the amebocyte lysate system may be used. Typical formulation additives may include, without limitation, about 100-300 mM NaCl, about 10-100 mM divalent cations (e.g.,  $Mg^{2+}$  or  $Ca^{2+}$ ), biocompatible buffers, e.g., Tris (tris(hydroxy)aminomethane), to give a final pH of about 6.0 to about 8.0, and, if the lysate is to be freeze dried, then  
5 sugars, e.g., mannitol or dextran. It is contemplated that the choice of appropriate formulation additives may also be determined by routine experimentation.

Synthetic chromogenic substrates have been used to measure the level of endotoxin-activated pro-clotting enzyme in LAL prepared from the hemolymph of both *Tachypleus tridentatus* and *Limulus polyphemus* horseshoe crabs (Iwanaga *et al.* (1978)  
10 HEMOSTASIS 7: 183-188). During an LAL assay that uses a chromogenic substrate, the pro-clotting enzyme (a serine protease) in the LAL is activated by endotoxin and cleaves the substrate's peptide chain on the carboxyl side of arginine so as to release the chromogenic group from the substrate, thereby releasing a marker compound that can be easily detected by, for example, spectrophotometry. One advantage of using a synthetic  
15 chromogenic substrate in an LAL assay in place of a conventional LAL gelation test is that the amount of activated clotting enzyme can be quantified and correlated to endotoxin levels in the sample.

Any chromogenic substrate that is cleaved by the clotting enzyme of a hemocyte lysate may be used in the practice of the invention. U.S. Patent No. 5,310,657, for  
20 example, describes an exemplary chromogenic substrate having the formula  $R_1-A_1-A_2-A_3-A_4-B-R_2$ , where  $R_1$  represents hydrogen, a blocking aromatic hydrocarbon or an acyl group;  $A_1$  represents an L or D-amino acid selected from Ile, Val or Leu;  $A_2$  represents Glu or Asp;  $A_3$  represents Ala or Cys;  $A_4$  represents Arg; B represents a linkage selected from an ester and an amide; and  $R_2$  represents a chromogenic or fluorogenic group which  
25 is covalently attached to the C-carboxyl terminal of Arginine through the B linkage, the fluorogenic or chromogenic moiety being capable of being cleaved from the remainder of the chromogenic substrate to produce a chromogen or a fluorogen. An exemplary chromogenic substrate has the consensus sequence acetate-Ile-Glu-Ala-Arg-pNA, where pNA represents a para-nitroaniline group. U.S. Patent No. 4,188,264 describes a peptide

substrate with a structure consisting of L-amino acids in the sequence R<sub>1</sub>-Gly-Arg-R<sub>2</sub>, where R<sub>1</sub> represents an N-blocked amino acid and R<sub>2</sub> is a group that can be released by enzymatic hydrolysis to yield a colored compound, HR<sub>2</sub>. U.S. Patent No. 4,510,241 discloses a chromogenic peptide substrate, which differs from the previous substrate in  
5 that the Gly moiety is replaced in the sequence by Ala or Cys. Alternatively, the chromogenic substrate may contain a fluorophore, for example, 7-amino-4-methyl coumarin, 7-amino-4-trifluoromethyl coumarin, and 4-methoxy-2-naphthylamine.

Inhibition or enhancement of the assay occurs when substances in the test sample interfere with the hemocyte lysate reaction. Inhibition results in a longer reaction time,  
10 indicating lower levels of microbial contamination than may actually be present in the test sample. Enhancement results in shorter reaction time, indicating higher levels of microbial contamination than may actually be present in the test sample. To verify the lack of inhibition or enhancement, an aliquot of test sample (or a dilution of the test sample) is “spiked” with a known amount of an agent representative of the microbial  
15 contaminant to be measured. It is recommended that the microbial contaminant spike results in a final microbial contaminant concentration in the sample equal to the midpoint, on a log basis, between the microbial contaminant concentration of the highest and lowest standards in the standard curve. For example, in an assay with a standard curve spanning from 50 Endotoxin Units (EU)/mL to 0.005 EU/mL, samples should be spiked  
20 to contain a final microbial contaminant concentration of 0.5 EU/mL. In an assay with a standard curve spanning from 1 EU/mL to 0.01 EU/mL, the microbial contaminant spike should result in a final microbial contaminant concentration of 0.1 EU/mL.

The spiked sample is assayed in parallel with the unspiked sample. The resulting microbial contaminant concentration in the unspiked sample and the microbial  
25 contaminant recovered in the spiked sample then are calculated. The microbial contaminant recovered should equal the known concentration of the spike within about 25%. If the test sample (or dilution) is found to inhibit or enhance the reaction, the sample may require further dilution until the inhibition or enhancement is overcome. Initially, one may want to screen for inhibition or enhancement by testing 10-fold

dilutions of test sample. Once the approximate non-inhibitory or non-enhancing dilution is determined, the exact dilution can be found by testing two-fold dilutions around this dilution. The degree of inhibition or enhancement will be dependent upon the concentration of the test sample. If several concentrations of the same sample are to be  
5 assayed, it is necessary to establish performance characteristics for each concentration independently.

In fabricating the cartridge of the invention, it is helpful to combine the amebocyte lysate and chromogenic substrate with at least one resolubilizing agent, such as a sugar or salt, and at least one anti-flaking agent, such as a polymer, prior to drying the lysate onto  
10 the solid support.

The resolubilizing agent preferably stabilizes the lysate in the dried form and facilitates resolubilization of the reagents during the assay. Useful resolubilizing agents include, for example, mannitol, mannose, sorbitol, trehalose, maltose, dextrose, sucrose, and other monosaccharides and disaccharides. The hemocyte lysate and chromogenic  
15 substrate preferably contain from about 0.01% (w/v) to about 20% (w/v), more preferably from about 0.1% (w/v) to about 1.0% (w/v) of the resolubilizing agent prior to drying.

The anti-flaking agent is an agent that prevents or reduces the likelihood that the hemocyte lysate and/or chromogenic substrate becomes disassociated from a solid support in the form of a dry flake. The anti-flaking agent preferably also stabilizes the  
20 hemocyte lysate or chromogenic substrate in the dried form. Useful anti-flaking agents include, for example, one or more polymers, including, for example, polyethylene glycol, polyvinyl pyrrolidone, dextrans, mannitol, and proteins, for example, serum albumin. The lysate preferably contains from about 0.01% (w/v) to about 25% (w/v), more preferably from about 0.1% (w/v) to about 1.0% (w/v) of anti-flaking agent prior to drying.

In addition, it has been found that certain polymers reduce the formation of air bubbles (e.g., frothing) when the hemocyte lysate and/or chromogenic substrate are resolubilized. Useful anti-frothing agents include polyvinyl alcohol and polypropylene glycol. In order to reduce frothing, the lysate and/or chromogenic substrate may contain  
25

from about 0.01% (w/v) to about 10% (w/v), more preferably from about 0.1% (w/v) to about 1.0% (w/v) anti-frothing agent prior to drying.

An exemplary fabrication process for the cartridge is described with reference to Figure 5, in which Figure 5A represents a bottom half 2 of cartridge 1 and Figure 5B represents a top half 3 of cartridge 1. Once prepared, the two halves of the cartridge 1 are joined to one another by adhesive, solvent bonding, ultrasonic welding, snap fit joints, or the like.

In Figure 5A, the bottom half 2 of the cartridge 1 defines one half of each conduit 8' (each having a first region 14' and a second region 16'). During fabrication of the bottom half 2 of the cartridge 1, hemocyte lysate is applied to each first region 14' and chromogenic substrate is applied to each second region 16'. In Figure 5B, the top half 3 of the cartridge 1 defines one half of each conduit 8''. During fabrication of top half 3 of the cartridge 1, an agent representative of a microbial contaminant (i.e., a spike), for example, a preselected amount of an endotoxin, is applied to region 14''. Once the reagents have been applied to the respective top 3 and bottom 2 halves of the cartridge 1, the cartridge halves 2 and 3 then are dried under conditions that preserve the activity of the hemocyte lysate and permit reconstitution of the hemocyte lysate to produce active lysate. In order to preserve the activity of the reagents during drying, the cartridge halves 2 and 3 are placed in an environment having a temperature from about 4°C to about 40°C, more preferably, from about 10°C to about 35°C, more preferably, from about 15°C to about 30°C, and a relative humidity from about 0% to about 30%, more preferably, from about 2% to about 20%, more preferably from about 4% to about 10%. Preferred drying conditions include a temperature of about 25°C and a relative humidity of about 5%. An exemplary protocol for manufacturing a cartridge of the invention is provided in Example 1.

In an alternative approach, the hemocyte lysate may be dried via freeze drying under standard conditions, about -30°C to about -40°C under vacuum.



After drying, the two cartridge halves **2** and **3** are joined to one another to create an intact cartridge **1**. Figure **5C** is a cross-sectional view through Section A-A' in which the two halves of the conduit (namely **8'** and **8''**) together create an intact conduit **8**, wherein region **14'** of the bottom **8'** of each conduit contains immobilized hemocyte lysate **20** and region **14''** of the top **8''** of one conduit contains immobilized endotoxin **22**. Figure **5D** is a cross-sectional view through Section B-B' in which region **16'** of the bottom **8'** of each conduit contains immobilized chromogenic substrate **24**.

Figures **6A-6D** are illustrations of two exemplary cartridges **1** of the invention corresponding to Figures **6A-6B** and Figures **6C-6D**. In Figure **6A**, the cartridge **1** may have an alternative finger grip **5** as shown with the dashed line. Figures **6A** and **6B** illustrate that the optical cell **6** in the first cartridge **1** is substantially cylindrical in shape. In Figure **6C**, the cartridge **1** also has a similar finger grip **5** to that shown by the dashed line in Figure **6A**. Figures **6C** and **6D** illustrate that the optical cell **6** in second cartridge **1** is more elongate in shape. The elongate shape permits greater depth and rise of fluid for greater optical pathlength and proportionally greater detection sensitivity. In addition, it is contemplated that the top and bottom halves **2** and **3** of each cartridge **1** may comprise one or more male (female) members and one or more reciprocal and interfitting female (male) members to stack the unassembled cartridge halves one on top of the other, as well as provide mating alignment in the assembled state.

The dimensions of a particular cartridge **1** may vary depending upon the number and/or type of assays to be performed. However, in one embodiment, as shown schematically in Figure **6A**, for example, the cartridge **1** has a length of about 10.16 cm (4.00"), width of about 2.54 cm (1.00"), and a height of about 0.476 cm (.188"). The bore of the conduit **8** running from the fluid inlet port **4** to the optical cell **6** is about .127 cm (.050"), where the lysate is dried on a region **14** of the conduit **8** about 2.381 cm (.938") from the fluid inlet port **4**, and a chromogenic substrate is dried on a region **16** of the conduit **8** about 4.65 cm (1.831") from the fluid inlet port **4**. The optical cell **6** in this embodiment is dimensioned to accommodate about 25  $\mu$ L of sample.

### Specimen Collection and Preparation

The cartridge may be used to determine the level of microbial contamination in a fluid, for example, a fluid to be administered locally or systemically, for example, parenterally to a mammal, or a body fluid to be tested for infection, including, for  
5 example, blood, lymph, urine, serum, plasma, ascites fluid, lung aspirants, and the like. In addition, the cartridge may be used to determine the level or microbial contamination in a water supply, for example, a supply of drinking water. In addition, the cartridge may be used to determine the level of microbial contamination in a food product, pharmaceutical, or medical device.

10 In general, materials used to harvest, store, or otherwise contact a sample to be tested, as well as test reagents, should be free of microbial contamination, for example, should be pyrogen-free. Materials may be rendered pyrogen-free by, for example, heating at 250°C for 30 minutes. Appropriate precautions should be taken to protect depyrogenated materials from subsequent environmental contamination.

### 15 Representative Assays That Can Be Performed in the Cartridge

It is contemplated that a variety of hemocyte lysate assays may be used in the cartridge of the invention, such as, for example, the end point turbidometric assay, the kinetic turbidometric assay, the endpoint chromogenic assay, and the single-step kinetic assay. In addition, the cartridge of the invention may be used with the multi-step kinetic  
20 assay, as described herein.

#### 1. End Point Turbidometric Assay

The end point turbidometric assay is described in Prior (1990) *supra*, pp. 28-34. Briefly, the end point turbidimetric assay includes the steps of (i) solubilizing a hemocyte lysate with a sample to be analyzed, (ii) incubating the resulting mixture at a temperature  
25 of about 0° to about 40° C, preferably about 25° to about 40°C, for a predetermined time, and (iii) measuring the increase in turbidity as a result of coagulation, if any, using a conventional coagulometer, nephrometer, or spectrophotometer.

Referring to Figure 2A, in order to perform an endpoint turbidometric assay in a cartridge 1, a sample is moved, for example, to a first region 14 of the conduit 8 containing the hemocyte lysate, where it is solubilized, for example, by cycling between forward and reverse pump action. The sample-lysate mixture then is moved to optical cell 6 for measurement of an optical property, for example, turbidity, using an optical detector. Results from multiple assays, for example, two assays can be averaged. The optical density of the sample-lysate mixture at a certain predetermined time point may then be interpolated onto a predetermined standard curve, for example, showing turbidity values on the Y axis versus endotoxin concentration on the X axis, to give the concentration of the microbial contaminant in the sample.

## 2. Kinetic Turbidometric Assay

The kinetic turbidometric assay is described in Prior (1990) *supra*, pp. 28-34. Briefly, the kinetic turbidometric assay includes the steps of (i) solubilizing a hemocyte lysate with a sample to be analyzed, (ii) incubating the resulting mixture at a temperature of about 0° to about 40°C, preferably about 25° to about 40°C, over a predetermined time range, and (iii) measuring a time required for either a turbidity change caused by coagulation to reach a pre-selected value or a ratio in change of the turbidity, using a conventional coagulometer, nephelometer, or spectrophotometer.

Referring to Figure 2A, in order to perform a kinetic turbidometric assay in a cartridge 1, a sample is, for example, moved to a first region 14 of the conduit 8 containing the hemocyte lysate, where it is solubilized, for example, by cycling between forward and reverse pump action. The sample-lysate mixture then is moved to optical cell 6 for measurement of an optical property, for example, turbidity, by measuring the absorbance or transmittance properties of the sample-lysate mixture using an optical detector. The detector may determine how long it takes for each optical property to exhibit, for example, a 5% drop in optical transmittance. Results from multiple assays, for example, two assays can be averaged. The resulting values may then be interpolated onto a predetermined standard curve, for example, showing time for a preselected change

in transmittance on the Y axis versus endotoxin concentration on the X axis, to give the concentration of the contaminant in the sample.

### 3. Endpoint Chromogenic Assay

The endpoint chromogenic assay is described in Prior (1990) *supra*, pp. 28-34,  
5 and U.S. Patent Nos. 4,301,245 and 4,717,658. Briefly, the endpoint chromogenic assay includes the steps of (i) solubilizing a hemocyte lysate preparation with a sample to be analyzed, (ii) incubating the resulting mixture at a temperature of about 0°C to about 40°C, preferably about 25°C to about 40°C, for a predetermined time, (iii) contacting a test device containing chromogenic substrate with the incubated sample-lysate mixture,  
10 (iv) adding a reaction inhibitor, and (v) measuring, e.g., by colorimetric change, a substance released from the synthetic substrate by protease activity.

Referring to Figure 2A, in order to perform an endpoint chromogenic assay in a cartridge 1, a sample is moved, for example, to a first region 14 of the conduit 8 containing the hemocyte lysate, where it is solubilized, for example, by cycling between  
15 forward and reverse pump action. Following a predetermined incubation period, the sample-lysate mixture then is moved, for example, by pump action to a second region 16 of the conduit 8 containing the chromogenic substrate, where it is solubilized, for example, by cycling between forward and reverse pump action. The sample-lysate-substrate mixture then is moved to a third region containing a reaction inhibitor.  
20 Afterwards, the sample-lysate-substrate mixture then is moved to optical cell 6 for measurement of an optical property, for example, the absorbance or transmittance properties of the sample by an optical detector. The optical property of the sample-lysate-substrate mixture at a certain predetermined time point may then be interpolated onto a predetermined standard curve, for example, showing absorbance, optical density, or  
25 transmittance on the Y axis versus endotoxin concentration on the X axis, to give the concentration of the microbial contaminant in the sample.

4. Single-Step Kinetic Assay

A single-step kinetic assay, for example, a single step-chromogenic assay, is described in U.S. Patent No. 5,310,657. Briefly, the kinetic chromogenic assay includes the steps of (i) simultaneously solubilizing a hemocyte lysate with a sample to be  
5 analyzed and a chromogenic substrate, (ii) incubating the resulting mixture at a temperature of about 0° to about 40°C, preferably about 25° to about 40°C, over a predetermined time range, and (iii) measuring a time required for a colorimetric change to reach a pre-selected value or a ratio in change of the colorimetric readout, using a conventional spectrophotometer.

10 Referring to Figure 2A, in order to perform a kinetic chromogenic assay in a cartridge 1, a sample is moved, for example, by pump action, to a first region 14 of the conduit 8 containing both the hemocyte lysate and chromogenic substrate, where it is solubilized, for example, by cycling between forward and reverse pump action. The sample-lysate-substrate mixture then is moved to optical cell 6 for measurement of an  
15 optical property for example, the absorbance or transmittance properties of the sample by an optical detector. The detector may determine how long it takes for each optical property to exhibit, for example, a 5% drop in optical transmittance. Results from multiple assays, for example, two assays can be averaged. The resulting values may then be interpolated onto a predetermined standard curve, for example, showing the time for a  
20 preselected change in absorbance or transmittance (as the case may be) on the Y axis versus endotoxin concentration on the X axis, to give the concentration of the contaminant in the sample.

Of the above methods, the endpoint chromogenic assay and the single-step kinetic chromogenic assays currently are considered the most rapid, sensitive, and economic  
25 assays for the detection of microbial contaminants, for example, endotoxin. However, both assays have their limitations. The endpoint chromogenic assay is rapid (about 15 minutes) but typically can only detect endotoxin concentrations in a range that is limited to about one log range (for example, about 1 EU/mL to about 0.1 EU/mL), with a sensitivity of about 0.1 EU/mL. Although the single-step kinetic chromogenic assay

measures endotoxin concentrations in a wider range of about 5 logs (for example, about 5 to about 0.05 EU/mL) with a high sensitivity of about 0.05 EU/mL, this method can be quite slow to run (about 30 minutes). Furthermore, neither the endpoint chromogenic assay nor the single-step kinetic chromogenic assay are readily adaptable for in-field performance. The multi-step kinetic assay overcomes the limitations in the endpoint chromogenic assay and the kinetic chromogenic assay.

#### 6. Multi-step Kinetic Assay

As will be discussed in more detail, the cartridge may also be used to perform a multi-step kinetic assay. The various steps involved in the multi-step kinetic assay are shown schematically in Figure 7. The assay is initiated by combining the sample to be tested with a volume of a hemocyte lysate to produce a sample-lysate mixture. The mixture then is incubated for a predetermined period of time. The sample-lysate mixture then is contacted with a substrate, for example, a chromogenic substrate, to produce a sample-lysate-substrate mixture. Thereafter, the time in which a preselected change in an optical property (for example, a specific change in an absorbance value or a specific change in a transmission value) is measured. The presence and/or amount of microbial contaminant may be then determined by interpolating the measured time against a pre-calibrated standard curve, for example, a standard curve showing the time to make a preselected change in optical property (absorbance or transmittance) on the Y axis versus endotoxin concentration on the X axis.

The standard curve may be created, for example, by adding increasing amounts of an agent, for example, an endotoxin, glucan, or other microbial cell wall component, in a blank sample, for example, pyrogen-free water. The time for which a preselected change in an optical property, for example, a preselected increase in absorbance or a preselected decrease in transmittance, is determined for each concentration of the microbial cell wall component. The various time measurements to achieve a standard change in optical property then are plotted as a function of the microbial cell wall component concentration. In general, the concentration of endotoxin is inversely proportional to the time necessary to achieve the standard change in optical property. The standard curve can

then be used to assess the presence and/or amount of microbial contaminant in the sample of interest. The calculation of standard curves is provided in Example 2.

As will be apparent to one skilled in the art, the relative amounts of hemocyte lysate and substrate can be adjusted to ensure that effective amounts of these two components are present in the sample-lysate-substrate mixture at the end of the assay. The final amount of hemocyte lysate protein in the assay is from about 1  $\mu\text{g}$  to about 500  $\mu\text{g}$ , preferably about 50  $\mu\text{g}$ . The final amount of the substrate, for example, the chromogenic substrate in the assay is from about 1  $\mu\text{g}$  to about 50  $\mu\text{g}$ , preferably about 6.5  $\mu\text{g}$ . The determination of the concentration and composition of the substrate, for example, the chromogenic substrate, is considered to be within the level of skill in the art.

The final volume of the resulting sample-lysate-substrate mixture can be based on the requirements of the optical detector used to measure the change in optical property of the sample. The ratio of volumes between the sample, lysate, and substrate can be readily established by those of ordinary skill in the art. Depending on the relative volumes of the sample, lysate, and substrate in the sample-lysate-substrate mixture, the concentration of the other components of the assay can be adjusted to maintain the final concentrations in the operable range, as described herein.

Referring to Figure 2A, to perform the multi-step kinetic assay in a cartridge 1 of the invention, a sample is first moved, for example, by pump action, to a first region 14 containing the hemocyte lysate, where it is mixed and incubated for a predetermined period of time. The sample-lysate mixture then is moved, for example, by pump action, to the second region 16 containing the substrate, for example, the chromogenic substrate, where it is solubilized. The sample-lysate-substrate mixture then is moved to optical cell 6, for a measurement of an optical property. The time intervals required for mixing and incubating steps are preprogrammed for optimal sensitivity and microbial contaminant concentration range.

Figures 8A and 8B show an exemplary cartridge and hand held optical detector useful in the practice of the invention. Figure 8A shows the cartridge being introduced

into the detector and Figure 8B shows the cartridge inserted into the detector with the fluid inlet ports still exposed to the user.

Figure 9 is a graph showing the changes in optical properties that can be generated in a cartridge using a multi-step kinetic assay. The dashed line represents changes in transmittance of the sample over time. The solid line represents changes in absorbance of the sample over time. Referring to Figures 2A and 9, optical properties were monitored after an aliquot of sample was added to the fluid inlet port 4 at a time of 0 seconds. After 60 seconds, the sample was drawn to region 14, where it was mixed with hemocyte lysate disposed on region 14 for 60 seconds (represented by vertical zig zag lines from T=60 to T=120). The amplitude of the mixing (length of the vertical lines) is determined such that the sample is repeatedly moved over the entire region 14. Thereafter, the sample-lysate mixture was incubated from 120 to 480 seconds at region 14 without further mixing. After 480 seconds, the sample was drawn to a chromogenic substrate at region 16 and the sample-lysate mixture combined with the chromogenic substrate from the period shown from 480 to 540 seconds (represented by vertical zig zag lines from T=480 to T=540). The amplitude of the mixing (length of the vertical lines) is determined such that the sample is repeatedly moved over the entire region 16. The resulting sample-lysate-substrate mixture then was drawn to optical cell 6 and the optical properties (absorbance and transmittance) measured for the period from 540 second to 1440 seconds.

Using the initial absorbance or transmittance readings of the mixture, the time required for the absorbance or transmittance to change by an arbitrary amount (Reaction Time) is determined. The amount of microbial contaminant in the sample then is determined by comparing the Reaction Time for the sample against a predetermined standard curve.

A spiked sample is assayed in parallel with the unspiked sample. The microbial contaminant concentration in the unspiked sample and the microbial contaminant recovered in the spiked sample can be compared to determine the presence of interference, such as an inhibitor or an enhancer of the reaction, as previously described.



Although the multi-step assay may be performed in a cartridge of the type discussed above, it may be employed in a variety of other formats, for example, within the well of a microtiter plate. Exemplary assays performed in the well of a microtiter plate are discussed in Example 3. Multiple samples of various test fluids, as well as  
5 spiked samples and the series of control samples making up a standard curve, may be placed in the wells of the microplate. Fixed amounts of hemocyte lysate and then substrate are added to each of the wells, preferably using an automated system, such as a robot, and the plate processed by a microplate reader, which can be programmed to sequentially read the absorbance of each well in a repetitive fashion.

10 In addition, it is contemplated that the cartridges, glucan-specific assays, and the multi-step kinetic assays can be used to detect the presence and/or amount of a ligand of interest in a test sample. For example, by adapting the assay format as appropriate, it is possible to detect the presence and/or amount of any ligand of interest, for example, a drug, toxin, protein, metabolite, in a sample. An exemplary ligand assay performed in the  
15 well of a microtiter plate is discussed in Example 7. By way of example, a binder for a ligand of interest, for example, an antibody or antigen binding fragment thereof, is immobilized on the surface of a solid support, for example, in a cartridge or a microplate. The binder is then pre-loaded or pre-bound with a complex comprising the ligand of interest coupled or bound to lipopolysaccharide or glucan. Methods for conjugating  
20 glucan or lipopolysaccharide to a ligand are well known in the art. For example, Boutonnier *et al.* (2001) *INFECT. IMMUN.* 69:3488-3493, describe methods for conjugating lipopolysaccharide to tetanus toxoid (see also, Konadu *et al.* (1994) *INFECT. IMMUN.* 62:5048-5054; Kenne *et al.* (1982) *CARBOHYDR. RES.* 100:341-349). When a sample containing the ligand of interest is combined with the immobilized binder, the  
25 lipopolysaccharide or glucan-labeled ligand is displaced. The amount of displaced ligand can then be quantified by the extent of the lipopolysaccharide or glucan initiated reaction of a LAL preparation.

Using these principles, it is possible to create a cartridge for determining the presence and/or amount of a ligand of interest in a sample. In this type of format, the

binder for ligand is immobilized on a surface of the conduit downstream of the sample inlet port and upstream of the optical cell. The binder for ligand is then preloaded with a complex comprising the ligand of interest coupled or conjugated to, for example, lipopolysaccharide or glucan. A hemocyte lysate (for example, a Factor G-specific lysate for detecting displaced glucan, or a Factor C-specific lysate for detecting displaced lipopolysaccharide) is immobilized on a surface of the conduit downstream of the immobilized binder for ligand. When the sample of interest is applied to the sample inlet port, the sample passes the binder for ligand. To the extent that the sample contains the ligand, the ligand displaces the complex from the immobilized binder. The amount of displaced ligand can be measured by measuring a change in an optical property in the hemocyte lysate. This change in optical property can then be correlated with the amount of the ligand of interest in the sample.

### **EXAMPLES**

Practice of the invention will be more fully understood from the following examples, which are presented herein for illustrative purposes only, and should not be construed as limiting the invention in any way.

#### **Example 1: Preparation of the Chromogenic Assay Cartridge**

An exemplary cartridge shown in Figures 2 was prepared as follows. Referring to Figure 5A, the LAL and chromogenic substrates were applied to regions 14' and 16', respectively, of conduit 8' of the bottom half 2 of the cartridge 1 using a Hamilton Microlab 540B Dispenser (Hamilton Company, Reno, NV). Briefly, 4-5.0  $\mu$ L of 10 mg/mL Endosafe LAL (Charles River Endosafe, Charleston, SC) containing 1% mannitol (Osmitrol, Baxter, Deerfield, IL) and 0.1% dextran (MW 10-100K, Sigma-Aldrich, St. Louis, MO), was applied to regions 14'. 4-5.0  $\mu$ L of (1.3 mg/mL) chromogenic substrate Ile-Glu-Ala-Arg-pNA Chromogenix S-2423 (Instrumentation Laboratories, Milan, Italy) containing 1 % polyvinyl alcohol (PVA) (MW 7K-30K, Sigma-Aldrich, St. Louis, MO), was applied to regions 16'. The bottom half 2 of the cartridge 1 was dried under a

controlled temperature of 25°C +/- 2°C and a humidity of 5% +/- 5% in a Linaire Environmental Steady State & Stability Test Chamber (Linaire Environmental, Williamsport, PA) in a Puregas HF200 Heatless Dryer (MTI Puregas, Denver, CO) for 1 hour. Temperature and humidity was controlled by a Watlow Series 96 1/16 DIN  
5 Temperature Controller (Watlow Electric Manufacturing Company, St. Louis, MO).

Referring to Figure 5B, 5µl Endosafe CSE endotoxin (Charles River Endosafe, Charleston, SC) (“spike”) was applied to region 14” of the conduit 8” of the top half 3 of the cartridge 1. The top half 3 of the cartridge 1 was dried under a controlled temperature of 25°C +/- 2°C and a humidity of 5% +/- 5% for one hour, as described above.

10 Following fabrication, the two halves 2 and 3 were assembled such that regions 14’ and 14” were aligned one on top of the other, and the edges of the cartridge halves 2 and 3 ultrasonically sealed using a Dukane Model 210 Ultrasonic Sealer (Dukane Corporation, St. Charles, IL) under the control of a Dukane Dynamic Process Controller (Dukane Corporation, St. Charles, IL).

15 The resulting cartridge 1 was labeled to identify the lot number of the cartridge, in order to later identify the standard curves used to quantify the microbial contaminant in the sample. The sealed, labeled cartridge 1 then was placed into a laminated foil pouch along with a desiccant such as silica gel or molecular sieve. The foil pouch was purged with nitrogen gas and then sealed with a PAC Model PV-G Vacuum Impulse Heat Sealer  
20 (Packaging Aids Corporation, San Rafael, CA).

#### Example 2: Cartridge-based Assays

This example demonstrates that a cartridge of the invention can be used to measure the amount of a microbial contaminant by an endpoint chromogenic assay, kinetic chromogenic assay and a multi-step chromogenic assay.

25 Figure 8 shows how the cartridge of the invention may be used with a portable hand-held optical detector. Figure 8A shows the cartridge in the process of being inserted into the optical detector. Figure 8B shows the cartridge inserted fully into the optical

detector, however, the fluid inlet ports of the cartridge are still accessible to the user. This configuration permits the user to apply one or more samples of interest to the exposed fluid inlet ports even though the optical cells of the cartridge are located in place within the optical detector.

5 (I) Endpoint Chromogenic Assay in a Cartridge

A cartridge was prepared essentially as described in Example 1, with the exception that no spikes were added to the conduits. Samples of endotoxin standards of 1.0 EU/mL, 0.5 EU/mL, 0.25 EU/mL, and 0.125 EU/mL were prepared and 25  $\mu$ L of each were pipetted into one of four cartridge fluid inlet ports. The portable optical  
10 detector maintained a temperature of 37°C for the duration of the test. The portable optical detector was programmed to perform a series of steps. The portable optical detector first drew each endotoxin sample into the region of a conduit that contained dried LAL, and mixed and incubate the sample with the LAL for 120 seconds (T=60 to T=180). The portable optical detector then drew the endotoxin sample-LAL mixture to  
15 the region in the conduit containing dried chromogenic substrate, and mixed the sample-LAL mixture with the substrate for 5 seconds. The portable optical detector then drew the sample-LAL-substrate mixture to the optical cell. The portable optical detector then recorded absorbance data from each of the four channels. After about 10 minutes (about 780 seconds), the test was ended by reading the last absorbance (optical density) value.  
20 The absorbance value curves for each endotoxin sample are shown in Figure 10A. The absorbance values generated at 780 seconds were recorded and plotted as a function of endotoxin concentration (Figure 10B) to give a standard curve. This standard curve can be used to determine the concentration of endotoxin in a sample of interest, when the sample is treated in the same manner as the standards.

25 (II) Single-Step Kinetic Chromogenic Assay in a Cartridge

A cartridge was prepared essentially as described in Example 1, with the exceptions that the chromogenic substrate was applied to the same region as the LAL and that no spikes were added to the conduits. Samples of endotoxin standards of 5.0

EU/mL, 0.5 EU/mL, and 0.05 EU/mL were prepared, and 25  $\mu$ L of the 5 EU/mL standard was pipetted into two fluid inlet ports of the cartridge. The portable optical detector maintained a temperature of 37°C for the duration of the test. The portable optical detector was programmed to perform a series of steps. The portable optical detector first  
5 drew each endotoxin sample into the region of a conduit that contained both dried LAL and chromogenic substrate, and mixed and incubated the sample with the LAL and substrate for 30 seconds. The portable optical detector then drew the sample-LAL-substrate mixture to the optical cell. The portable optical detector began recording absorbance for both samples. The assay was repeated for the 0.5 EU/mL and 0.05  
10 EU/mL standards.

The plots of recorded data for each endotoxin standard are shown in Figures 11A, 11B, and 11C. The time taken for the optical density of each endotoxin standard-LAL-substrate mixture to reach an optical density of 0.05 was recorded as the onset time for each standard. A plot of the log of the endotoxin concentration (X axis) vs. the log of the  
15 onset times (Y axis) provides a kinetic standard curve (Figure 11D). This standard curve can be used to determine the concentration of endotoxin in a sample of interest when the sample is treated in the same manner as the standards.

### (III) Multi-step Kinetic Assay in a Cartridge

A cartridge was prepared essentially as described in Example 1, with the  
20 exception that no spikes were added to the conduits. Samples of endotoxin standards of 1.0 EU/mL, 0.5 EU/mL, 0.25 EU/mL, and 0.125 EU/mL were prepared and 25  $\mu$ L of each were pipetted into one of four cartridge fluid inlet ports. The portable optical detector maintained a temperature of 37°C for the duration of the test. The portable detector was programmed to perform a series of steps. The portable optical detector first  
25 drew each endotoxin sample into the region of a conduit that contained dried LAL, and mixed and incubate the sample with the LAL for 240 seconds. The portable optical detector then drew the endotoxin sample-LAL mixture to the region in the conduit containing dried chromogenic substrate and mixed the sample-LAL mixture with the

substrate for 20 seconds. The portable optical detector then drew the sample-LAL-substrate mixture to the optical cell. The portable optical detector began recording absorbance data from each of the four channels. The time taken for the optical density of each endotoxin standard-LAL-substrate mixture to reach an optical density of 0.05 was recorded as the onset time for each standard (see, Table 1). A plot of the log of the endotoxin concentrations (X axis) versus the log of the onset times (Y axis) provides a kinetic standard curve (Figure 12). This standard curve can be used to determine the concentration of endotoxin in a sample when the sample is treated in the same manner as the standards.

10

TABLE 1

Endotoxin Concentration, EU/mL	5.0	0.5	0.05	0
Onset Time or Reaction Time (seconds)	30	166	222	300

Example 3: Microplate-based Assays

This example demonstrates that a multi-well microtiter plate can be fabricated so that the amount of a microbial contaminant (bacterial endotoxin in this example) can be measured by an endpoint chromogenic assay, a single-step kinetic assay and a multi-step assay.

(I) Single-step Kinetic Assay on a Microplate

A single-step kinetic chromogenic assay was performed as follows. Briefly, 50  $\mu$ L of a control microbial contaminant of interest (e.g., 5 EU/mL, 0.5 EU/mL, or 0.05 EU/mL of Endosafe Control Standard Endotoxin (CSE), Charles River Endosafe, Charleston, SC), was added to one or more wells of a 96 well plate. 50  $\mu$ L of 5 mg/mL Endosafe LAL (Charles River Endosafe, Charleston, SC) and 5  $\mu$ L of 1.3 mg/mL Chromogenix S-2423 chromogenic substrate (Instrumentation Laboratories, Milan, Italy) were added to each well and incubated at 37°C for 32 minutes. The optical density of the mixture in each well was monitored by a spectrophotometer, such as a Sunrise micro

25

plate reader (Tecan, Research Triangle Park, NC). The time taken for each standard to change 0.1 absorbance units was determined (“onset time”). The results are summarized below in Table 2. A plot of the log of the endotoxin concentrations (X axis) versus the log of the onset times (Y axis) provides a kinetic standard curve (Figure 13). This standard curve may be used to measure the concentration of endotoxin in a sample of interest, when the sample is treated in the same manner as the standards.

TABLE 2

Standard	Concentration, EU/mL	Time to Onset OD/ Max Reaction Time (seconds)	Mean Time to Reach Onset OD (seconds)	Standard Deviation	(RT) CV%	Calculated Value
STD1	5.0	229.4 / 221.4	225.4	4.0 / 4.0	1.77	>4.9853
STD2	0.5	472.2 / 465.0	468.6	3.59 / 3.59	0.77	0.5029
STD3	0.05	982.5 / 977.1	979.8	2.68 / 2.68	0.27	<0.0499
CTRL1	0	>990.0 / >990.0	>990.0	0.00 / 0.00	0.00	<0.0500

(II) Endpoint Assay on a Microplate

10 An endpoint chromogenic assay using the reagents in Example 3(I) was performed as follows. Briefly, 50 µL of a control microbial contaminant of interest (e.g., endotoxin) at a concentration of e.g., 5 EU/mL, 0.5 EU/mL, or 0.05 EU/mL, was added to one or more wells of a 96-well plate. 50 µL of hemocyte lysate (5 mg/mL) was added to one or more wells containing the standard and incubated at 37°C for 5 minutes. 5 µL of chromogenic substrate (1.3 mg/mL) was added to one or more wells and incubated at 15 37°C for 5 minutes. The reaction then was stopped by the addition of 100µL of 50% acetic acid to each well. The optical density of the mixture in each well was measured at 405 nm by a spectrophotometer, such as a Sunrise micro plate reader (Tecan, Research Triangle Park, NC). The absorbance of each sample at 780 seconds (shown in Table 3) was recorded. The absorbance of each sample (Y axis) was plotted versus the endotoxin concentration (X axis) to provide a standard curve (shown in Figure 14). This standard 20

curve may be used to measure the concentration of endotoxin in a sample of interest, when the sample is treated in the same manner as the standards.

TABLE 3

Standard	Concentration, EU/mL	Optical Density (OD)	Mean Optical Density	Standard Deviation	(OD) CV%	Calculated Value
STD1	1.2	0.7505 / 0.8105	0.7805	0.03 / 0.03	3.84	1.1876
STD2	0.6	0.4135 / 0.3825	0.3980	0.02 / 0.02	3.89	0.6319
STD3	0.3	0.1655 / 0.1615	0.1635	0.00 / 0.00	1.22	0.2912
STD4	0.15	0.0655 / 0.0525	0.0590	0.01 / 0.01	11.02	0.1393
CTRL1	0	0.0115 / -0.0115	0.0000	0.01 / 0.01	0.00	<0.1500

5 (III) Multi-step Kinetic Assay on a Microplate

A multi-step kinetic chromogenic assay using the same reagents as in Example 3(I) was performed as follows. Briefly, 50 µL of a control microbial contaminant of interest, for example, endotoxin, at the following concentrations of 5 EU/mL, 0.5 EU/mL, or 0.05 EU/mL, was added to one or more wells of a 96-well plate. 50 µL of hemocyte lysate (5 mg/mL) was added one or more wells and incubated at 37°C for 3 minutes. 5 µL of chromogenic substrate (1.3 mg/mL) was added to one or more wells containing the standard and incubated at 37°C for 16.5 minutes. The optical density of the mixture in each well was monitored by a spectrophotometer, such as a Sunrise micro plate reader (Tecan, Research Triangle Park, NC). The time taken for each microbial contaminant standard and/or sample to change 0.1 absorbance units was determined (“onset time”). The results are summarized in Table 4. A plot of the log of the endotoxin concentrations (X axis) versus the log of the onset times (Y axis) provides a kinetic standard curve (Figure 15). This standard curve may be used to measure the concentration of endotoxin in a sample of interest, when the sample is treated in the same manner as the standards.



TABLE 4

Standard	Concentration, EU/mL	Time to Onset OD/ Max Reaction Time (seconds)	Mean Time to Reach Onset OD (seconds)	Standard Deviation	(RT) CV%	Calculated Value
STD1	5.0	172.1 / 166.4	169.3	2.84 / 2.84	1.68	>5.1505
STD2	0.5	393.8 / 387.7	390.7	3.03 / 3.03	0.78	0.4712
STD3	0.05	851.7 / 843.3	847.5	4.17 / 4.17	0.49	<0.0515
CTRL1	0	>990.0 / >990.0	>990.0	0.00 / 0.00	0.00	<0.0500

**Example 4: Preparation and Testing of Glucan-Specific LAL**

5 As shown in Figure 1, crude LAL reacts with both endotoxin (lipopolysaccharide) and glucan, so that cell wall material from both Gram negative bacteria and yeast/mold cells activate the coagulation cascade. In the case of Gram negative bacteria, coagulation is mediated through the Factor C cascade. In the case of yeast and mold, coagulation is mediated through the Factor G cascade. In order to determine which of the two

10 contaminants is present in a sample, or what proportion of each might comprise a mixed contamination, LAL was produced under conditions to render it specific for glucan, as described below.

Amebocyte blood was mixed 6:1 with 0.05% Tween 20 (Sigma, St. Louis, MO), 150 mM NaCl and centrifuged in a Sorval model RC-3B centrifuge at 3,000 rpm for 5

15 minutes. The pelleted cells were washed with 0.01% Tween 20, 150 mM NaCl and centrifuged at 3,000 rpm for 5 minutes. The pelleted, washed cells were divided into three aliquots, resuspended in either lipopolysaccharide-free water, lipopolysaccharide-free 1 M NaCl, or lipopolysaccharide-free 2 M NaCl. The cells in each pellet were lysed by sonication for 1-2 minutes. Cell debris was removed by centrifugation at 4,000 rpm

20 for 10 minutes, and the resulting supernatant was harvested and used directly in coagulation experiments.

Figure 16 provides a graphical representation of kinetic coagulation reactions using various concentrations of lipopolysaccharide (Figure 16A) or glucan (Figure 16B)

in a microtiter plate multi-step kinetic assay, such as that described in Example 3 (III) with either standard LAL (row A of each figure) or glucan-specific LAL (rows B and C in each figure).

Figure 16A illustrates the reactivity of the lysates with lipopolysaccharide (Charles River Endosafe) serially diluted 1:10 from column 1 to column 5. The lipopolysaccharide concentrations in columns 1 through 5 are: 10 ng/ml, 1 ng/ml, 100 pg/ml, 10 pg/ml, and 0, respectively.

Figure 16B illustrates the reactivity of the lysates with glucan (Charles River Endosafe, Charleston, SC) serially diluted 1:10 from column 1 to column 5. Glucan concentrations from column 1 through 5 are: 100 µg/ml, 10 µg/ml, 1 µg/ml, 100 ng/ml, and 0, respectively.

In each figure, row A shows the reactivity with standard LAL prepared from cells lysed in pyrogen-free water (i.e., reactive with both glucan and lipopolysaccharide). In each figure, row B is glucan-specific LAL produced by lysing cells in 1 M NaCl. In each figure, row C is glucan-specific LAL produced by lysing cells in 2M NaCl. Each graph represents the change in optical density or absorbance, as shown on the Y axis, over time, as shown on the X axis. For example, in Figure 16A, the graph shown in row A, column 1, shows the change in absorption over time when 10 ng/ml lipopolysaccharide is added to standard LAL.

The results demonstrate that when lipopolysaccharide is added to standard lysate, the lipopolysaccharide activates the lysate to produce an increase in absorbance (see, Figure 16A, row A). However, the rate of absorbance change decreases as less lipopolysaccharide is added to each sample. The results demonstrate, however, that there is substantially no change in absorbance as lipopolysaccharide is added to each glucan-specific lysate (see, Figure 16A, rows A and B).

The results also demonstrate that when glucan is added to standard lysate, the glucan (like the lipopolysaccharide) activates the lysate to produce an increase in

absorbance (see, Figure 16B, row A). However, the rate of absorbance change decreases as less glucan is added to each sample. In contrast to the situation when lipopolysaccharide was added, the glucan activates each glucan-specific lysate to produce an increase in absorbance over time (see, Figure 16B, rows B and C).

5            These results demonstrate that it is possible to produce a glucan-specific lysate using the protocol described herein.

Example 5:    Testing of Glucan-Specific LAL in a Cartridge-based Multi-Step Assay

10            Glucan-specific LAL was prepared by lysing amebocytes in 2M NaCl as described in Example 4 and tested in a cartridge-based multi-step kinetic assay, such as that described in Example 2(III).

15            Samples containing glucan at 100 µg/ml, 10 µg/ml, 1 µg/ml, 100 ng/ml, and 0, were incubated with the glucan-specific lysate for 4 minutes. The resulting mixture was mixed with a chromogenic substrate, acetate-Ile-Glu-Ala-Arg-pNA and the change in absorbance at 405 nm was measured over time. The time required to reach an onset optical density of 0.05 was collected (see, Table 5) and the log of the glucan concentration was plotted versus the log of the time to reach onset O.D. (see, Figure 17).

TABLE 5

Glucan Standard (µg/ml)	Onset Time (seconds)	Calculated µg/ml glucan activity
100	207	80.5
10	400	15.4
1	1300	0.80
0.1	> 1800	< 0.36
Negative control	> 1800	< 0.36

20

The data show that it is possible to produce a standard curve of glucan concentration using a cartridge-based multi-step assay. The standard curve can then be used to determine the concentration of glucan in a sample of interest, when the sample is treated in the same manner as the standards.

25

Example 6: Means of Reading an LAL Reaction Using Fluorescent Substrates

A cartridge-based multi-step kinetic LAL assay, such as that described in Example 2(III), was performed using the fluorogenic substrate: Glu-Gly-Arg-AMC (Enzyme Systems Products, Livermore, CA). The cartridge was modified to include a long-pass filter placed between the sample and the light sensor such that light having the excitation wavelength (390 nm) was blocked but yet the emissions at a wavelength of 460 +/- 25 nm were able to pass through and be detected by the sensor.

Using the device, fluorescence data (expressed as Relative Fluorescence Units) were collected from samples containing 10 EU (Endotoxin units), 1 EU, and 0.1 EU. The changes in Relative Fluorescence Units over time are presented in Figure 18. The results demonstrate that the multi-step kinetic assays can measure the amount of endotoxin in sample using a fluorescent substrate.

Example 7: Measurement of Lipopolysaccharide-labeled Ligand Using an Immobilized Antibody

The microtiter plate multi-step kinetic LAL assay, such as that described in Example 3(III) can also be used to measure the concentration of a ligand. Briefly, an antibody that binds a ligand of interest is immobilized onto the surface of a well of a microtiter plate. Then, the ligand binding sites of the antibody are preloaded with a complex comprising the ligand coupled to lipopolysaccharide. When a sample containing the ligand is exposed to the immobilized antibody, the lipopolysaccharide-labeled ligand is displaced and quantified by reaction with LAL. In this example, the ligand detected was fluorescein. An anti-fluorescein antibody was immobilized onto the surface of a well and was pre-loaded with fluorescein-labeled lipopolysaccharide. When samples containing fluorescein-labeled antibody were exposed to the immobilized antibody, the fluorescein-lipopolysaccharide conjugate was released from the immobilized antibody and measured by reactivity with LAL.

Briefly, rabbit anti-fluorescein antibody (Virostat, Portland ME) was diluted 1/4000 in CAPS buffer, pH 10.2 (Sigma, St. Louis, MO). 25 µl of the diluted antibody

was added to the wells of a high binding plate (Corning 25801) and incubated for 1 hour at 37 °C. The plate was washed 4 x 100 µl per well with TTBS (0.1% Tween, 100 mM Tris buffered saline) using a multipipettor. The wells were blocked by adding 150 µl per well of gelatin diluent and stored at 4 °C overnight. The wells then were washed with 3 x 100 µl per well with TTBS. Lipopolysaccharide (LPS)-fluorescein (FITC) conjugate (List Biological Laboratory, Campbell, CA) was diluted from 1 µg/ml to 10 pg/ml using 10-fold dilutions in 0.1 M Tris. 0.1 M Tris was used as a control. 75 µl of diluted LPS/FITC conjugate was added per well and incubated for 30 minutes at 37°C. The wells then were washed with 3 x 100 µl per well with 0.1 M Tris.

Fluorescein-labeled goat anti-chicken antibody (Southern Biotech. Associates Inc 6100-02, Birmingham, AL) was diluted 1/50, 1/150, and 1/450 in 0.1 M Tris. 0.1 M Tris was used as a control. 75 µl of fluorescein-labeled goat anti-chicken antibody was added per well and incubated for 30 minutes at 37 °C. 50 µl was removed from each well and transferred to a clean 96-well plate (Falcon, 353072, Becton Dickenson, Franklin Lakes, NJ). 50 µl Endochrome K (Charles River Endosafe), which is a 1:1 mixture of LAL substrate and LAL lysate, was added to each well. The plate was read at time intervals (e.g., kinetically) at 405 nm at 37 °C for 60-90 minutes (Min OD: 0, Max OD: 0.8, Onset OD: 0.1).

To find the proper concentration of LPS/FITC to be pre-bound to the immobilized antibody, increasing concentrations were applied to the plate from 10 pg/mL to 1 µg/mL. Then dilutions of the fluorescein-labeled antibody ligand were exposed to the immobilized antibody and the displaced LPS/FITC measured as EU equivalents. As shown in Figure 19, at the 1 µg/mL level of LPS/FITC, an EU proportional to the ligand was achieved (far left). Accordingly, by using this type of format, it is possible to detect the presence and/or measure the amount of a ligand of interest in a test sample.

### Equivalents

The invention may be embodied in other specific forms without departing from the spirit or essential characteristics thereof. The foregoing embodiments are therefore to

be considered in all respects illustrative rather than limiting on the invention described herein. Scope of the invention is thus indicated by the appended claims rather than by the foregoing description, and all changes that come within the meaning and range of equivalency of the claims are intended to be embraced therein.

5

Incorporation by Reference

All publications and patent documents cited in this application are incorporated by reference in their entirety for all purposes to the same extent as if the entire contents of each individual publication or patent document was incorporated herein.